

Brochure

Overlays



Blake
Medical Distribution

1-866-662-3459
www.blakemedical.ca



Sleep Comfortably & Reduce Pressure Injuries

The Geo-Matrix™ Gel overlay is constructed of 1.5" thick dry polymer gel. The overlay is designed for use over an existing foam mattress. This overlay is indicated for prevention and maintenance of skin health and the treatment of early stage, non-complex pressure ulcers. Each overlay comes with a four-way stretch, fluid proof, fire retardant (Cal.Bulletin 117) anti-microbial cover.


The Geo-Matrix™ Gel columns support a certain amount of weight. When this threshold has been met, the columns collapse to allow for gentle immersion. Peak pressure points are largely reduced through the columns ability to load share. As bony prominences contact the columns, the load is distributed through their neighbouring columns.

Multiple gel columns create segmentation. This segmentation greatly reduces shear by decreasing the interface pressure when users move or transfer across the surface.

Product Code: GMO(+size), Example GMO3680

Standard Sizes: 36x80, 39x80, 42x80, 48x80, 54x80 and 60x80

Warranty: 2 years for Gel Overlay, 1 year for Cover

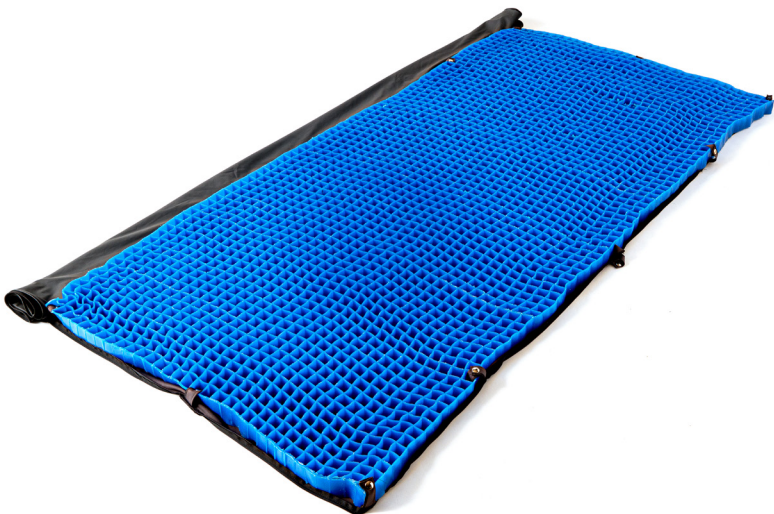

The Foam overlay is constructed of 2" thick gel infused memory foam. The overlay is designed for use over an existing foam mattress. This overlay is indicated for prevention and maintenance of skin health and the treatment of early stage, non-complex pressure ulcers. Each overlay comes with a four-way stretch, fluid proof, fire retardant (Cal.Bulletin 117) anti-microbial cover.

Product Code: GIO(+size), Example GIO3680

Standard Sizes: 36x80, 39x80, 42x80, 48x80, 54x80 and 60x80

Warranty: 2 years for Foam Overlay, 1 year for Cover

