




Proudly made
and manufactured
in Canada!

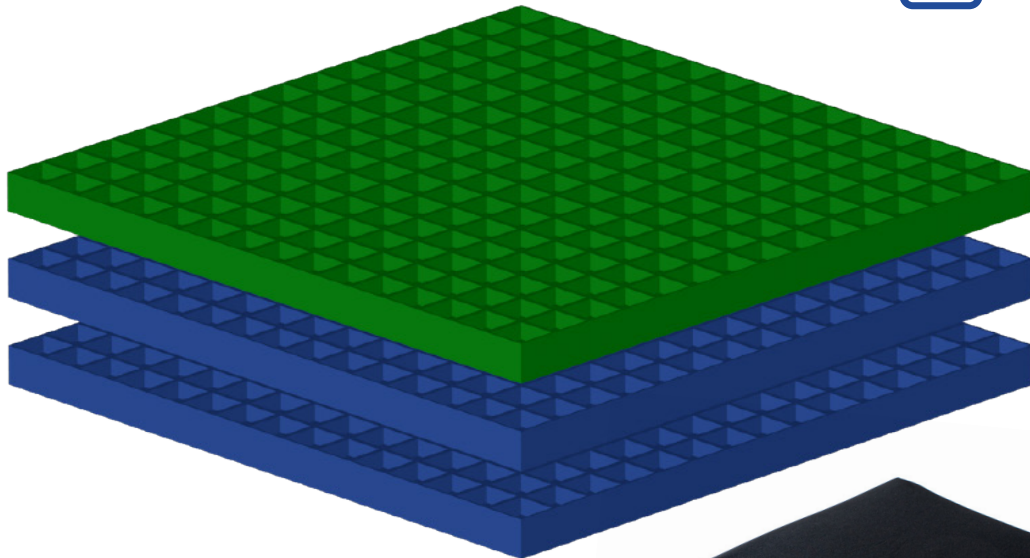

Product Specifications

CONSTRUCTION

Two layers of 2.5cm high Blue gel and one layer of 2.5cm high Green gel.

COVERS

Two outer cover options: Breathable Spacetex fabric with incontinent protection OR Elastatex soft touch outer incontinent cover.

Inner cover: Elastatex breathable, wipeable, antibacterial, fluid proof cover.

STANDARD SIZES

35x35, 35x38, 35x40, 35x43, 35x48, 35x50, 38x35, 38x38, 38x40, 38x43, 38x45, 38x48, 38x50, 40x35, 40x38, 40x40, 40x43, 40x45, 40x48, 40x50, 40x53, 40x55, 43x35, 43x38, 43x40, 43x43, 43x45, 43x48, 43x50, 43x55, 45x30, 45x33, 54x35, 45x38, 45x40, 45x43, 45x45, 45x48, 45x50, 45x53, 45x55, 48x35, 48x38, 48x40, 48x43, 48x45, 48x48, 48x50, 48x53, 48x55, 50x35, 50x38, 50x40, 50x43, 50x45, 50x48, 50x50, 50x53, 50x55, 53x35, 53x38, 53x40, 53x43, 53x45, 53x48, 53x50, 53x53, 53x55, 55x35, 55x38, 55x40, 55x43, 55x45, 55x48, 55x50, 55x53, 55x55, 58x40, 58x43, 58x45, 58x48, 58x50, 58x53, 58x55, 60x40, 60x43, 60x45, 60x48, 60x50, 60x53, 60x55, 60x55

CUSTOM SIZES

All dimensions available

PRODUCT CODE

GG"WxD" eg. GG4545

WEIGHT CAPACITY


Recommended maximum 115kg

The **Geo-Matrix™ Gold** cushion is a three layer all gel cushion. It is commonly used with power or tilt chairs; often used in long-term care where frequent cleaning is required. These individual gel layers can be machine laundered and / or disinfected easily.

The Gold cushion provides maximum pressure redistribution for people with high risk of skin breakdown and/or existing wounds, limited mobility and/or postural control.

The Geo-Matrix™ Gold is ideal for users whom require vibration dampening along with extreme skin protection. The Geo-Matrix Gold is extremely low maintenance, as it requires no daily calibration and has zero risk of deflation or leaks.

Blake Medical esteemed distributor:

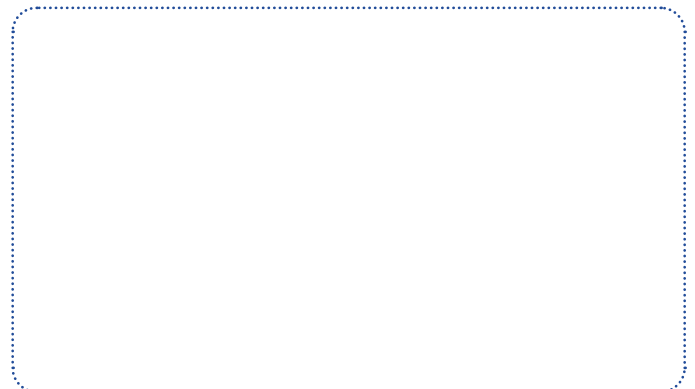

GEL CUSHION
Geo-Matrix Gold



WASHING INSTRUCTIONS


Covers: 60°C in washing machine. Hang to air dry. Do not bleach.

Geo-Matrix™ gel: Submerge and/or wipe down gel using any cleaner that does not contain bleach and allow to air dry.

OR - Add individual gel sheets to a front loading washing machine for a full cycle.

For more specific cleaning instructions please visit our website or user's manual.

All of Blake Medical products are free from phthalates and latex. This includes all medical fabrics, Geo-Matrix™ gel, moulded polyurethane, standard polyurethane and specialty foam. All cushions meet EN-1021 part 1 & 2 requirements for fire retardancy.

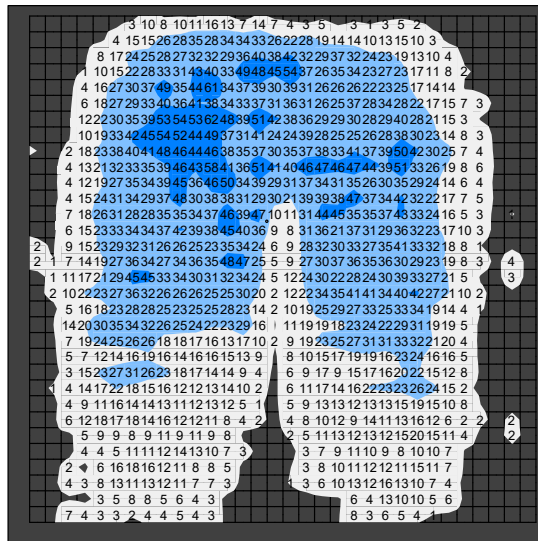

Client: Male 85kg

ADDITIONAL BENEFITS

- Cushion remains 2-3 degrees cooler under seated body temperature to assist with microclimate management.
- Minimizes shear through column segmentation.
- Geo-Matrix™ Gel does not change structure in climates of extreme heat nor extreme cold.
- Geo-Matrix™ gel is recyclable through Blake Medical.

GEL CUSHION
Geo-Matrix Gold

Advantages of Geo-Matrix™ Gel

DURABLE	COMFORTABLE	DISINFECT-ABLE	FRICTION & SHEAR REDUCING	INCREASES SIT TIMES
MICROCLIMATE MANAGEMENT	PRESSURE REDISTRIBUTING	VIBRATION DAMPENING	WASHABLE	WON'T PUNCTURE

