




Proudly made
and manufactured
in Canada!

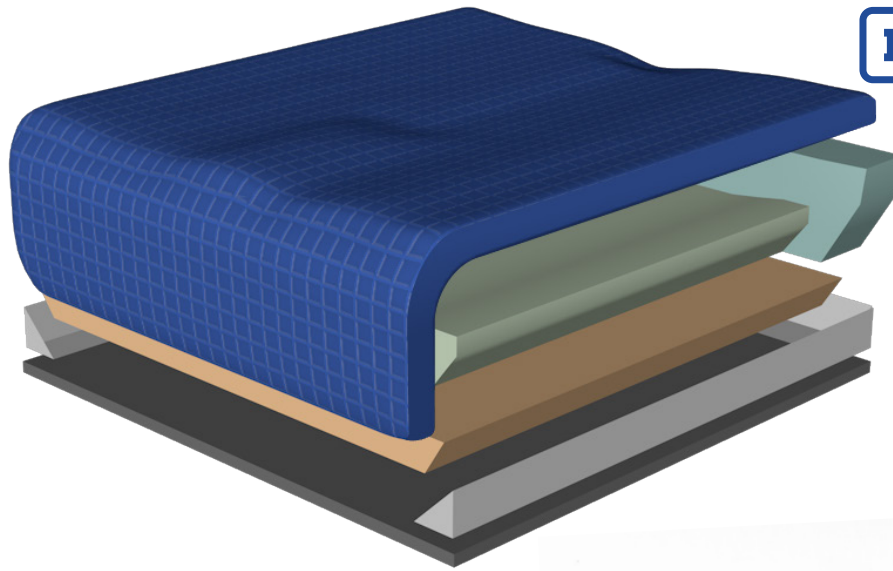

Product Specifications

CONSTRUCTION

Closed cell foam base, dual high-density foam core with a firm sacral shelf and cascading Pillow Top gel topper

COVERS

Two outer cover options: Breathable Spacetex fabric with incontinent protection OR Elastatex soft touch outer incontinent cover.

Inner cover: Elastatex breathable, wipeable, antibacterial, fluid proof cover.

STANDARD SIZES

35x35, 35x38, 35x40, 35x43, 35x48, 35x50, 38x35, 38x38, 38x40, 38x43, 38x45, 38x48, 38x50, 40x35, 40x38, 40x40, 40x43, 40x45, 40x48, 40x50, 40x53, 40x55, 43x35, 43x38, 43x40, 43x43, 43x45, 43x48, 43x50, 43x55, 45x30, 45x33, 54x35, 45x38, 45x40, 45x43, 45x45, 45x48, 45x50, 45x53, 45x55, 48x35, 48x38, 48x40, 48x43, 48x45, 48x48, 48x50, 48x53, 48x55, 50x35, 50x38, 50x40, 50x43, 50x45, 50x48, 50x50, 50x53, 50x55, 53x35, 53x38, 53x40, 53x43, 53x45, 53x48, 53x50, 53x53, 53x55, 55x35, 55x38, 55x40, 55x43, 55x45, 55x48, 55x50, 55x53, 55x55, 58x40, 58x43, 58x45, 58x48, 58x50, 58x53, 58x55, 60x40, 60x43, 60x45, 60x48, 60x50, 60x53, 60x55, 60x55

CUSTOM SIZES

All dimensions available

PRODUCT CODE

GH"WxD" eg. GH4545


WEIGHT CAPACITY

Recommended maximum 125kg

The Geo-Matrix™ Hybrid cushion is often used to provide comfort for users who require gentle positioning and contouring with excellent pressure distribution and vibration control.

The cascading gel top layer increases comfort and reduces pressure behind the knees. The Geo-Matrix™ Hybrid is designed for moderate postural support.

Blake Medical esteemed distributor:

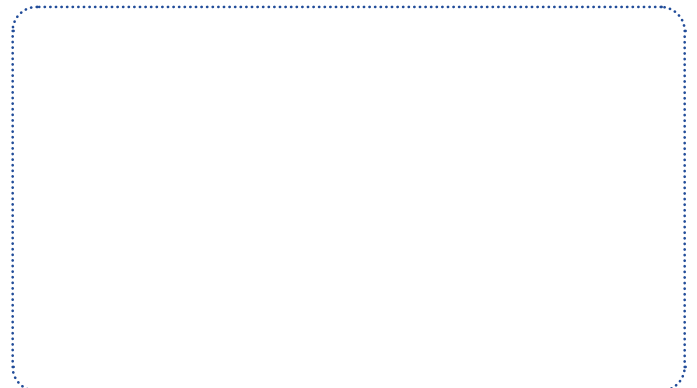

*Can be upgraded to 160kg weight capacity upon request

Blake Medical | 60 Prince Charles Road, Brantford ON, N3S 0H4 | blakemedical.ca | +1(905) 560-0773

GEL & FOAM CUSHION
Geo-Matrix Hybrid




WASHING INSTRUCTIONS

Covers: 60°C in washing machine. Hang to air dry. Do not bleach.

Geo-Matrix™ gel and foam: Wipe down gel using any cleaner that does not contain bleach and allow to air dry.

For more specific cleaning instructions please visit our website or user's manual.

All of Blake Medical products are free from phthalates and latex. This includes all medical fabrics, Geo-Matrix™ gel, moulded polyurethane, standard polyurethane and specialty foam. All cushions meet EN-1021 part 1 & 2 requirements for fire retardancy.

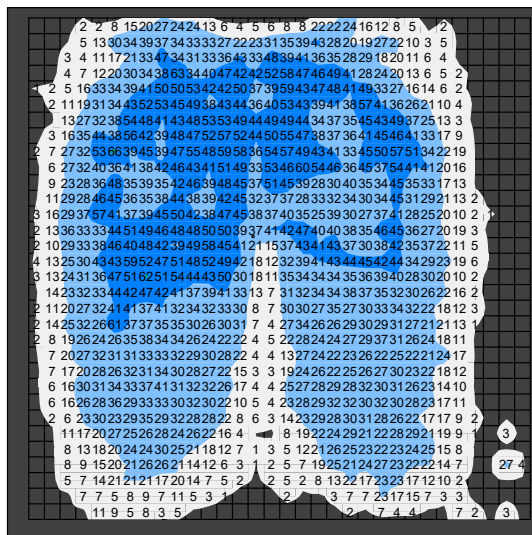

Client: Female 62kg

ADDITIONAL BENEFITS

- Cushion remains 2-3 degrees cooler under seated body temperature to assist with microclimate management.
- Minimizes shear through column segmentation.
- Geo-Matrix™ Gel does not change structure in climates of extreme heat nor extreme cold.
- Geo-Matrix™ gel is recyclable through Blake Medical.

Advantages of Geo-Matrix™ Gel

DURABLE	COMFORTABLE	DISINFECT-ABLE	FRICTION & SHEAR REDUCING	INCREASES SIT TIMES
MICROCLIMATE MANAGEMENT	PRESSURE REDISTRIBUTING	VIBRATION DAMPENING	WON'T PUNCTURE	

