




Proudly made  
and manufactured  
in Canada!

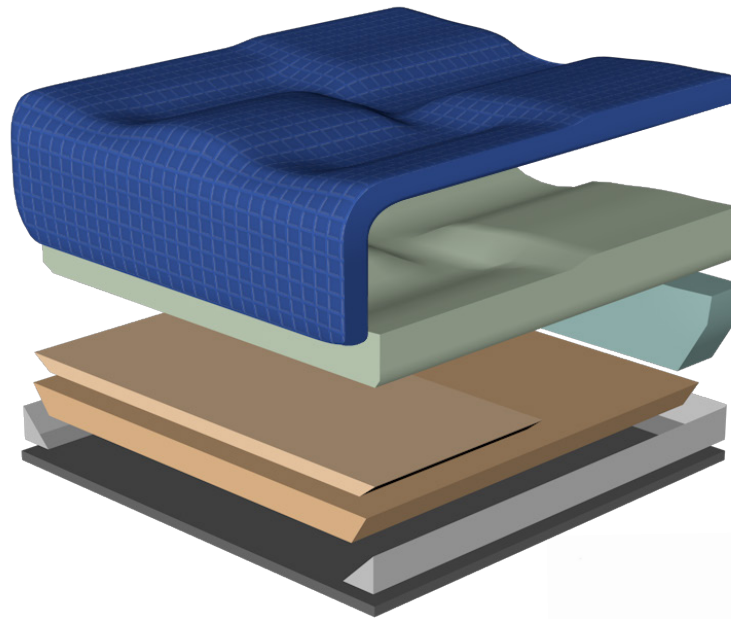

## Product Specifications

### CONSTRUCTION

Closed cell foam base, dual high-density foam core with a firm ischial and sacral shelf and cascading Pillow Top gel topper

### COVERS

Two outer cover options: Breathable Spacetex fabric with incontinent protection OR Elastatex soft touch outer incontinent cover.

Inner cover: Elastatex breathable, wipeable, antibacterial, fluid proof cover.

### STANDARD SIZES

35x35, 35x38, 35x40, 35x43, 35x48, 35x50, 38x35, 38x38, 38x40, 38x43, 38x45, 38x48, 38x50, 40x35, 40x38, 40x40, 40x43, 40x45, 40x48, 40x50, 40x53, 40x55, 43x35, 43x38, 43x40, 43x43, 43x45, 43x48, 43x50, 43x55, 45x30, 45x33, 54x35, 45x38, 45x40, 45x43, 45x45, 45x48, 45x50, 45x53, 45x55, 48x35, 48x38, 48x40, 48x43, 48x45, 48x48, 48x50, 48x53, 48x55, 50x35, 50x38, 50x40, 50x43, 50x45, 50x48, 50x50, 50x53, 50x55, 53x35, 53x38, 53x40, 53x43, 53x45, 53x48, 53x50, 53x53, 53x55, 55x35, 55x38, 55x40, 55x43, 55x45, 55x48, 55x50, 55x53, 55x55, 58x40, 58x43, 58x45, 58x48, 58x50, 58x53, 58x55, 60x40, 60x43, 60x45, 60x48, 60x50, 60x53, 60x55, 60x55, 60x55

### CUSTOM SIZES

All dimensions available

### PRODUCT CODE

GHD"WxD" eg. GHD4545


### WEIGHT CAPACITY

Recommended maximum 125kg

The **Geo-Matrix™ Hybrid Deep** cushion is often used to provide comfort for users who require optimal positioning and contouring with excellent pressure distribution and vibration control.

The cascading gel top layer increases comfort and reduces pressure behind the knees. The Geo-Matrix™ Hybrid Deep is designed for robust postural support.

Blake Medical esteemed distributor:

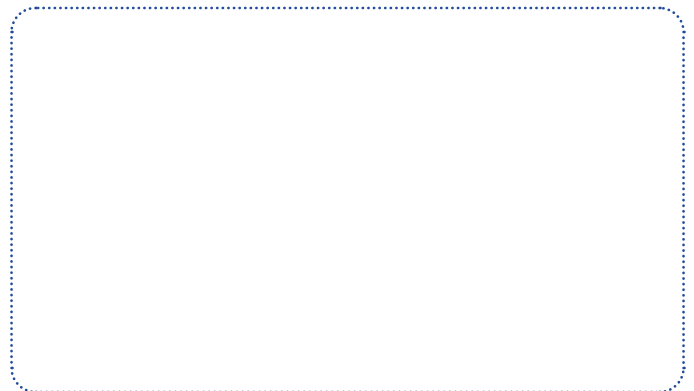

\*Can be upgraded to 160kg weight capacity upon request

Blake Medical | 60 Prince Charles Road, Brantford ON, N3S 0H4 | [blakemedical.ca](http://blakemedical.ca) | +1(905) 560-0773

# GEL & FOAM CUSHION Geo-Matrix Hybrid Deep




## WASHING INSTRUCTIONS

Covers: 60°C in washing machine. Hang to air dry. Do not bleach.

Geo-Matrix™ gel and foam: Wipe down gel using any cleaner that does not contain bleach and allow to air dry.

**For more specific cleaning instructions please visit our website or user's manual.**

All of Blake Medical products are free from phthalates and latex. This includes all medical fabrics, Geo-Matrix™ gel, moulded polyurethane, standard polyurethane and specialty foam. All cushions meet EN-1021 part 1 & 2 requirements for fire retardancy.

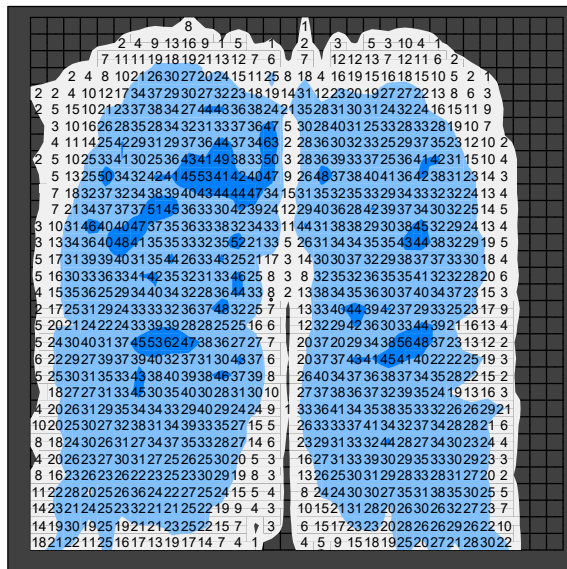

Client: Female 55kg

## ADDITIONAL BENEFITS

- Cushion remains 2-3 degrees cooler under seated body temperature to assist with microclimate management.
- Minimizes shear through column segmentation.
- Geo-Matrix™ Gel does not change structure in climates of extreme heat nor extreme cold.
- Geo-Matrix™ gel is recyclable through Blake Medical.

## Advantages of Geo-Matrix™ Gel

