




Proudly made
and manufactured
in Canada!

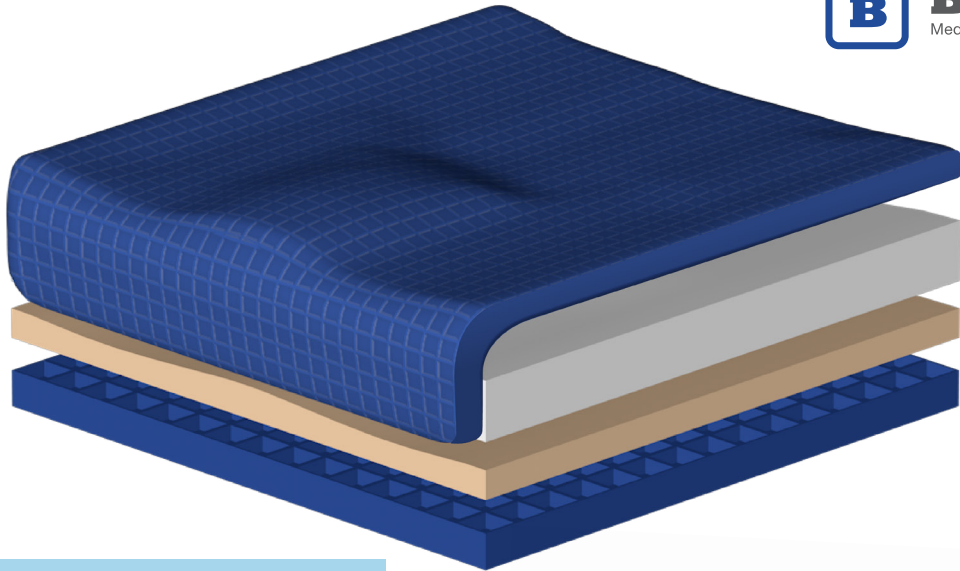

Product Specifications

CONSTRUCTION

A base layer of 2.5 cm thick Blue gel with a dual density contoured foam core and a top layer of 2 cm Pillow Top gel that cascades over the front of the cushion

COVERS

Two outer cover options: Breathable Spacetex fabric with incontinent protection OR Elastatex soft touch outer incontinent cover.

Inner cover: Elastatex breathable, wipe-able, antibacterial, fluid proof cover.

STANDARD SIZES

35x35, 35x38, 35x40, 35x43, 35x48, 35x50, 38x35, 38x38, 38x40, 38x43, 38x45, 38x48, 38x50, 40x35, 40x38, 40x40, 40x43, 40x45, 40x48, 40x50, 40x53, 40x55, 43x35, 43x38, 43x40, 43x43, 43x45, 43x48, 43x50, 43x55, 45x30, 45x33, 54x35, 45x38, 45x40, 45x43, 45x45, 45x48, 45x50, 45x53, 45x55, 48x35, 48x38, 48x40, 48x43, 48x45, 48x48, 48x50, 48x53, 48x55, 50x35, 50x38, 50x40, 50x43, 50x45, 50x48, 50x50, 50x53, 50x55, 53x35, 53x38, 53x40, 53x43, 53x45, 53x48, 53x50, 53x53, 53x55, 55x35, 55x38, 55x40, 55x43, 55x45, 55x48, 55x50, 55x53, 55x55, 58x40, 58x43, 58x45, 58x48, 58x50, 58x53, 58x55, 60x40, 60x43, 60x45, 60x48, 60x50, 60x53, 60x53, 60x55, 60x55

CUSTOM SIZES

All dimensions available

PRODUCT CODE

GP"WxD" eg. GP4545


WEIGHT CAPACITY

Recommended maximum 125 kg



The **Geo-Matrix™ Platinum** cushion is designed to provide comfort for people who require skin protection and gentle positioning with maximum pressure redistribution and increased stability. A base layer of gel provides vibration dampening, making this cushion perfect for powered wheelchairs. The Platinum is beneficial for those with existing skin issues or who are at higher risk for skin breakdown.

Blake Medical esteemed distributor:

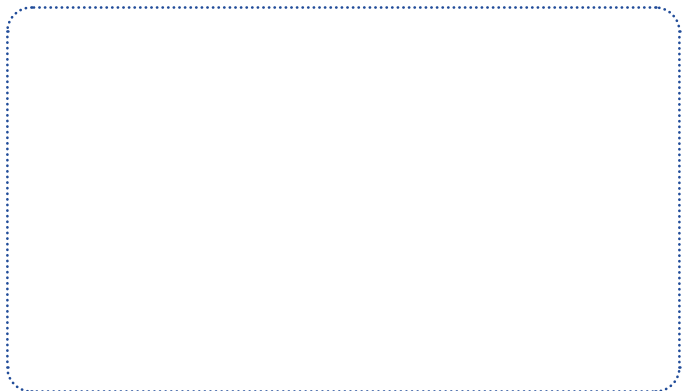

GEL & FOAM CUSHION Geo-Matrix Platinum




WASHING INSTRUCTIONS

Covers: 60°C in washing machine. Hang to air dry. Do not bleach.

Geo-Matrix™ gel and foam: Wipe down gel using any cleaner that does not contain bleach and allow to air dry.

For more specific cleaning instructions please visit our website or user's manual.

All of Blake Medical products are free from phthalates and latex. This includes all medical fabrics, Geo-Matrix™ gel, moulded polyurethane, standard polyurethane and specialty foam. All cushions meet EN-1021 part 1 & 2 requirements for fire retardancy.

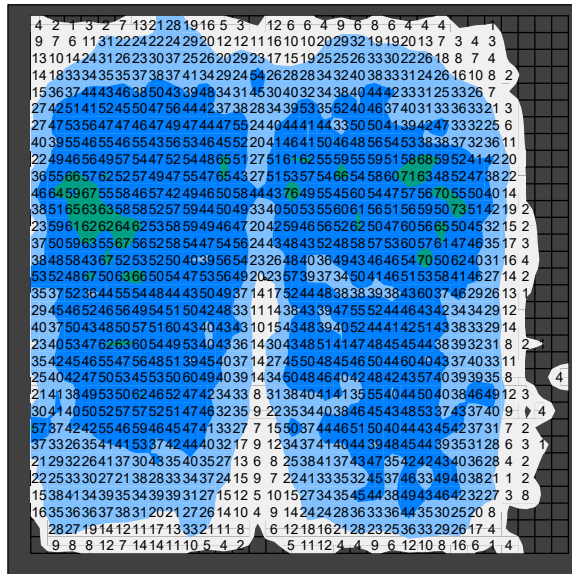

ADDITIONAL BENEFITS

- Cushion remains 2-3 degrees cooler under seated body temperature to assist with microclimate management.
- Minimizes shear through column segmentation.
- Geo-Matrix™ Gel does not change structure in climates of extreme heat nor extreme cold.
- Geo-Matrix™ gel is recyclable through Blake Medical.

Client: Male 102kg

Advantages of Geo-Matrix™ Gel

